



The World's Best Massage Chairs

Presidential Massage Chair

User 's Manual



*Make your life healthy
and beautiful every day!*



The World's Best Massage Chairs

Please read SAFETY INSTRUCTIONS carefully before using the chair.

INFINITY

Presidential Massage Chair

Function Introduction

1. Infinity Presidential features a new patented soft 3D massage mechanism. It is capable of reproducing kneading, knocking, shiatsu, and fast rubbing back massage like a live professional massage therapist.
2. The actual massage rollers are heated so they can perform heat massage on the neck, shoulder, and neck. Helps relax stressed muscles and promote blood circulation.
3. Super long rail design allows for multi-angle massage from neck to thigh. Amazing whole body massage!
4. Chair moves forward before massage, allowing you to put the chair close to the wall. Saves space and looks great!
5. One button puts the chair in full zero gravity or lie flat position, allowing you to relax your mind and body completely.
6. Ergonomic remote control uses a beautiful high definition LCD to show massage area, massage techniques, and system state in real time. Allows for easy selection of customized massage programs.
7. Chair is equipped with state of the art bluetooth technology allowing the option of operating the chair with phone or tablet and not using the remote at all.
8. Bluetooth wireless music player integrates music player technology with massage technology. This allows the massage to sync with your music allowing a personalized experience with every massage.
9. With omnibearing rollers, air pressure massage function, and foot heating system, your daily stress will melt away.
10. Armrests are equipped with pulse air pressure massage to promote blood circulation in the arms and relaxation.
11. Equipped with multiple sets of airbags in the waist and hip area.
12. Footrest can extend over 7 inches allowing for comfortable massage for different height people.

Bluetooth function: 1, download the Infinity APP from to the smart phone or tablet PC; 2, turn on the Bluetooth function at the smart phone or tablet PC 3, choose Infinity and connect it; 4, open the APP, then you can control the massage chair ,just like a remote control.

Bluetooth wireless music: After you connect smart phone or tablet PC with the chair, open the music player software, then you can play the music .

- ★ 1.Thank you for choosing Infinity Presidential Massage Chair.
- 2.Please read this manual carefully before using the massage chair.
- 3.Please keep this manual in a safe place for further reference.
- 4.Please pay attention to SAFETY INSTRUCTIONS.

★ Note: We reserve the right for design modification. It is subject to change without notice. The pictures in this manual are only for reference. Please refer to the actual product.

Content

Preparation	Safety instructions	1
	External structure parts	5
	Internal structure parts	6
	Mechanical components diagram	7
	Mechanical components list	8
Installation	Installation	9
Massage Operation	The adjustment before massage	13
	Manual control	15
	Operation instruction &function illustration	21
	After massage	24
Maintenance	Cleaning and maintenance	25
Others	Malfunction guide	26
	Technical data	27

Please keep this manual in a safe place for future reference.

Safety instructions

- Please read this manual carefully to make sure the installation and operation of the massage chair is correct.
- “Safety Instructions” is provided in both Warning and Caution for important safety information.

 Warning	This sign indicates warning of causing personal injury.
 Caution	This sign indicates possibility of damage to the chair.

Symbol meaning

Preparation

	This sign indicates forbidden behavior.		This sign indicates self disassembly is forbidden.
	This sign indicates must follow instruction to operate.		This sign indicates correct operation of cable.

- The following should be kept for any further references.

Safety instructions



Warning

- Osteoporosis patient.
- People with heart disease including embedded pacemaker or other medical electronic device inside the body.
- For pregnant women or people with health problems, please consult doctor before use.
- Patients with a high fever forbid to use this machine.
- Children should not use this chair.
- Please don't sleep in the chair.
- Do not use this chair when intoxicated.

People who are not suitable to use the chair

- Make sure there are no kids or pets or other objects around (under, behind or in front of) the chair at the moment of backrest, leg-rest up/down or leg-rest retraction/extension from reclining position.
- During operation, if you have body abnormal, please stop the operation and go to consult your physician.
- Please do not take strong massage on the head.
- Do not massage belly or knees, be careful when massage neck.
- Children or pets are not allowed to play on the chair especially during massage. Standing or sitting on the back or armrest is not allowed either.
- Those who move around difficulty or can't express themselves must be taken care of strictly and continuously during massage.
- Insert the plug of the chair into correct and suitable power socket. Make sure the plug is completely inserted into the socket to reduce short circuit or fire risk.
- Unplug the power when not in use, stop all massage functions and turn off on/off key before unplugging.

Accident/injury

- It is prohibited to use the massage chair while any items covered the massage chair ventilation holes.
- Massage cushions of the chair are not allowed to fold for storage.
- The machine is forbidden for no-load running for a long time.

Fire and electric shock

- Power cord, plug, socket can't be used in broken or loose contact condition.
- Don't use other electric power except rated voltage in case of fire, electric shock or any breakdown.
- Any broken, repaired, knotted, twisted power cord should be disposed of and replaced.

Accident

- If the chair works dysfunctional, or is broken, splashed with water or falls into water, please have it repaired by a professional repairman or send the chair to assigned service agent. Self-modification, disassembly or maintenance is forbidden.
- If power cord or plug is damaged, please don't replace it yourself, please contact us or appointed dealer for replacement.
- To avoid electric shock or bodily harm, please don't open or disassemble any part of the chair, there is not any attached part of the chair for maintenance by customers.
- Check to see if leather or cloth material on the back and other areas is broken or not. If broken, even though it is small, please stop using the chair, unplug the power and have it repaired

Disassembly and maintenance

Preparation

 Caution	
 <ul style="list-style-type: none"> Please operate the chair in a flat or level place Position the chair 4" from the wall to allow for sufficient space to recline up and down. The maximum suggested operating time is 20mins. 5mins for each position beginning with the slight-step massage. Massaging too long in one position is not recommended. Before sitting on the chair to avoid injury please inspect everything to confirm that the massage heads are in the correct position. Remote control should be put in the pocket when not in use. Do not put the remote control anywhere casually in case of remote control damage. Please make sure that the power switch is at the "OFF" position before inserting the plug into the socket. Turn off all controls before removing plug from outlet. Make sure the plug completely insert into the socket. Checking whether any damage on wire, make sure no water or damage on wire. Installing earth wire in correct way. The plug must be inserted into socket with earth ground.  	To avoid sudden falling. To make sure the back lying. Give your body time to recover. To avoid accidental injury. To avoid the accident when turn on the power. To reduce short circuits and any fire hazard. To avoid electric shock in any leakage or breakdown.
 <ul style="list-style-type: none"> Don't press with hands or step on the massage heads with feet. Don't put hands between the working massage heads. Don't put hands or feet in the gap of mechanical parts. Don't drop or embed anything into the slot of massage chair. Don't stand or sit or put anything on the backrest, armrest or leg-rest. Don't stand on the chair during operation. Don't move the chair during operation. Don't pull out the plug or turn off the power during operation. Don't put your hands or head between leg-rest or beneath. Don't use the massage chair while wearing hair ornaments. Make sure pockets are empty when massaging. Don't sit or press on the backrest when the chair is in reclining position. Don't use chair with bare skin or thin clothing. Can cause skin irritation. 	To avoid accidents or injury
 <ul style="list-style-type: none"> Keep the chair away from heat, wet, sharp things and corrosive materials. To use electric blanket or other heating products on the chair is forbidden. Don't put in the space of dust serious or corrosive gas environment, Don't put the chair in a small space places or use in the bad air circulation places. Keep the chair away from stove or other heating products. Don't put the chair directly under the sun or other high temperature places. Don't use the chair outdoors. 	To avoid fading and discoloration on upholstery or hardening on PU.

 Caution	
 <ul style="list-style-type: none"> Don't use or store the chair in the wet environment like nearing the swimming pool, or bathroom. Keep the chair and remote control dry all the time. 	To avoid leakage or electric shock.
 <ul style="list-style-type: none"> Make sure to install the earth wire in correct way (Incorrect install way may cause electric shock. If you don't know how to install, please find a professional electrician) Don't use plug adapter (If the enclosed plug doesn't match with socket, please don't change the plug by yourself and find a professional electrician to install a matched socket.) 	To avoid electric shock.
 <ul style="list-style-type: none"> Unplug the power before cleaning the chair. Unplug the power after using the chair. Unplug the power when any damage on the power. Hold the plug to unplug the power but not try to pull the wire to unplug the power. Please turn off the power switch and unplug the power after massage. 	To avoid electric shock or fire harzard.
 <ul style="list-style-type: none"> Unplug the power immediately when electricity goes off. 	To avoid accidents
	To avoid damage or accidents when power on again.

Grounding instruction


★ The chair is equipped with a cord having an grounding conductor. The plug must be plugged into an appropriate socket that is properly installed and matched with the plug.

Replace fuse


★ When the fuse is burned out, please turn off the switch and unplug the power. Replace the fuse with another of the same rating. Or you can find an electrician to replace the fuse.

 Earth wire is not allowed to be in these places	Gas pipe: Could cause fire. Phone wire and lightning rod: probably cause electric shock and fire. Water pipe: Plastic will not provide correct ground.
--	--


External structure

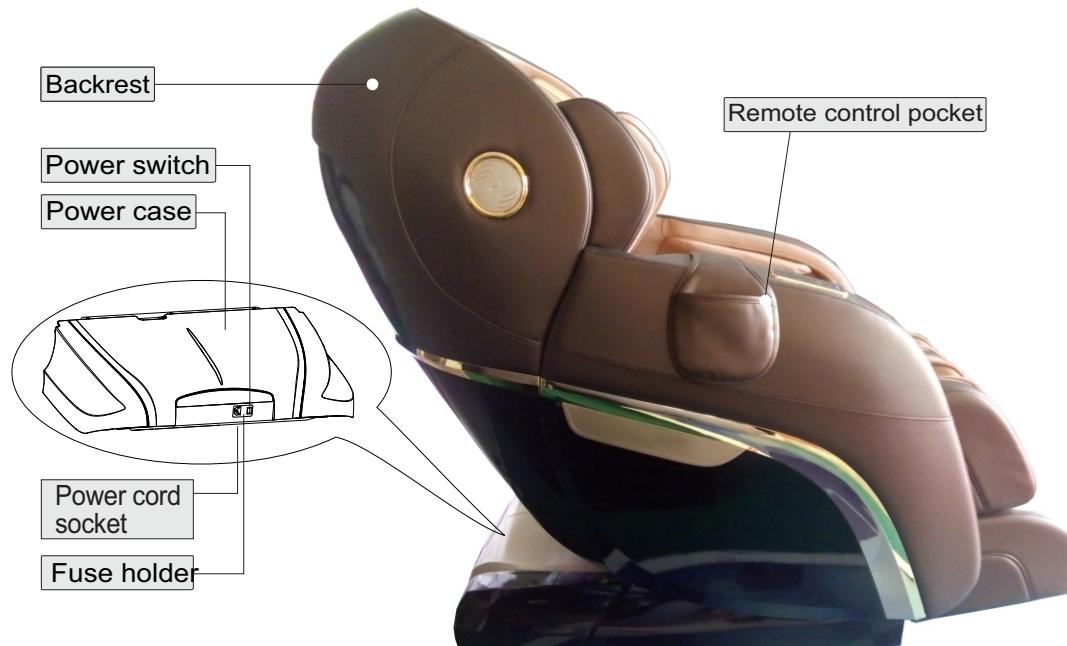

Internal structure

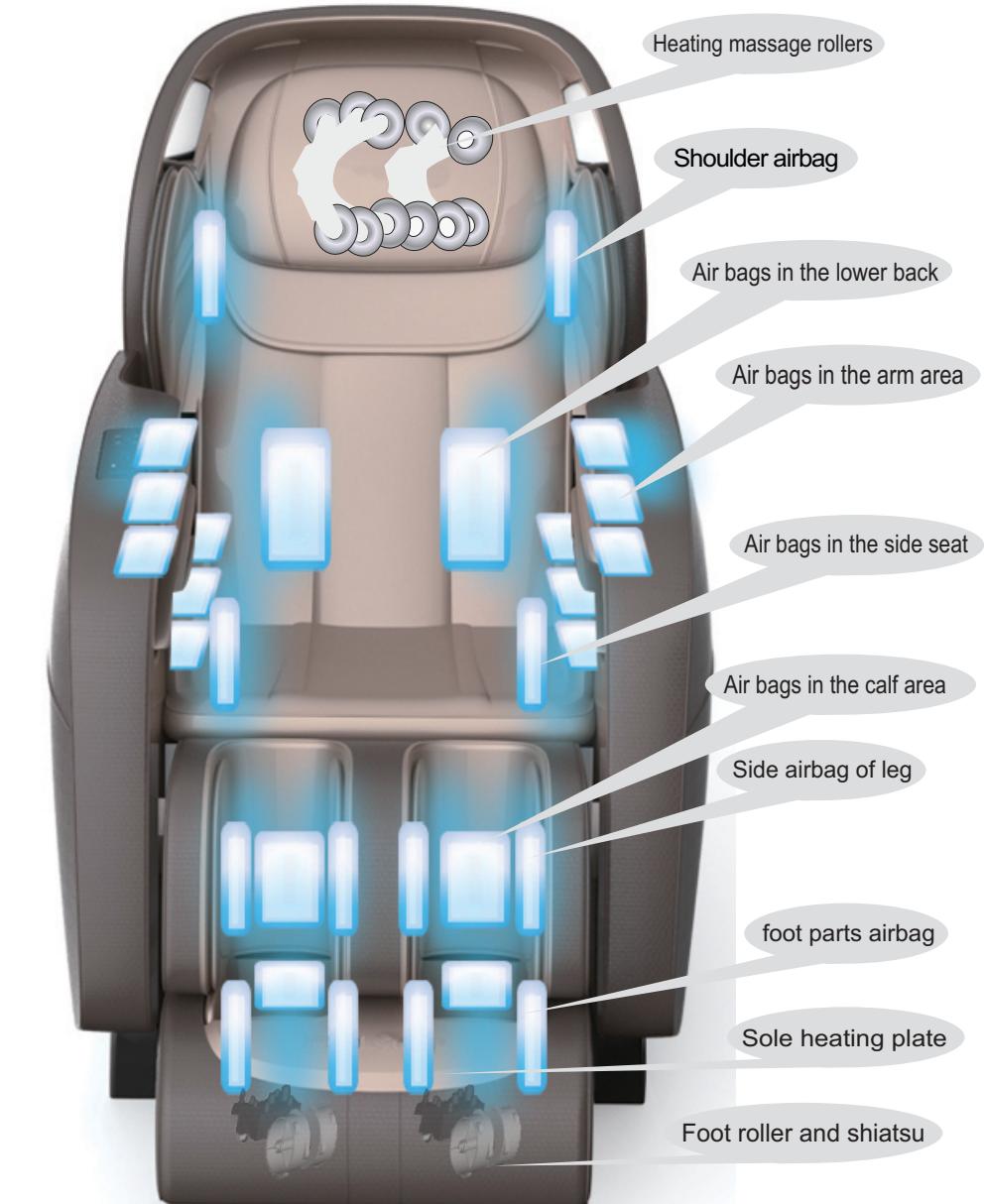

Explosion Diagram

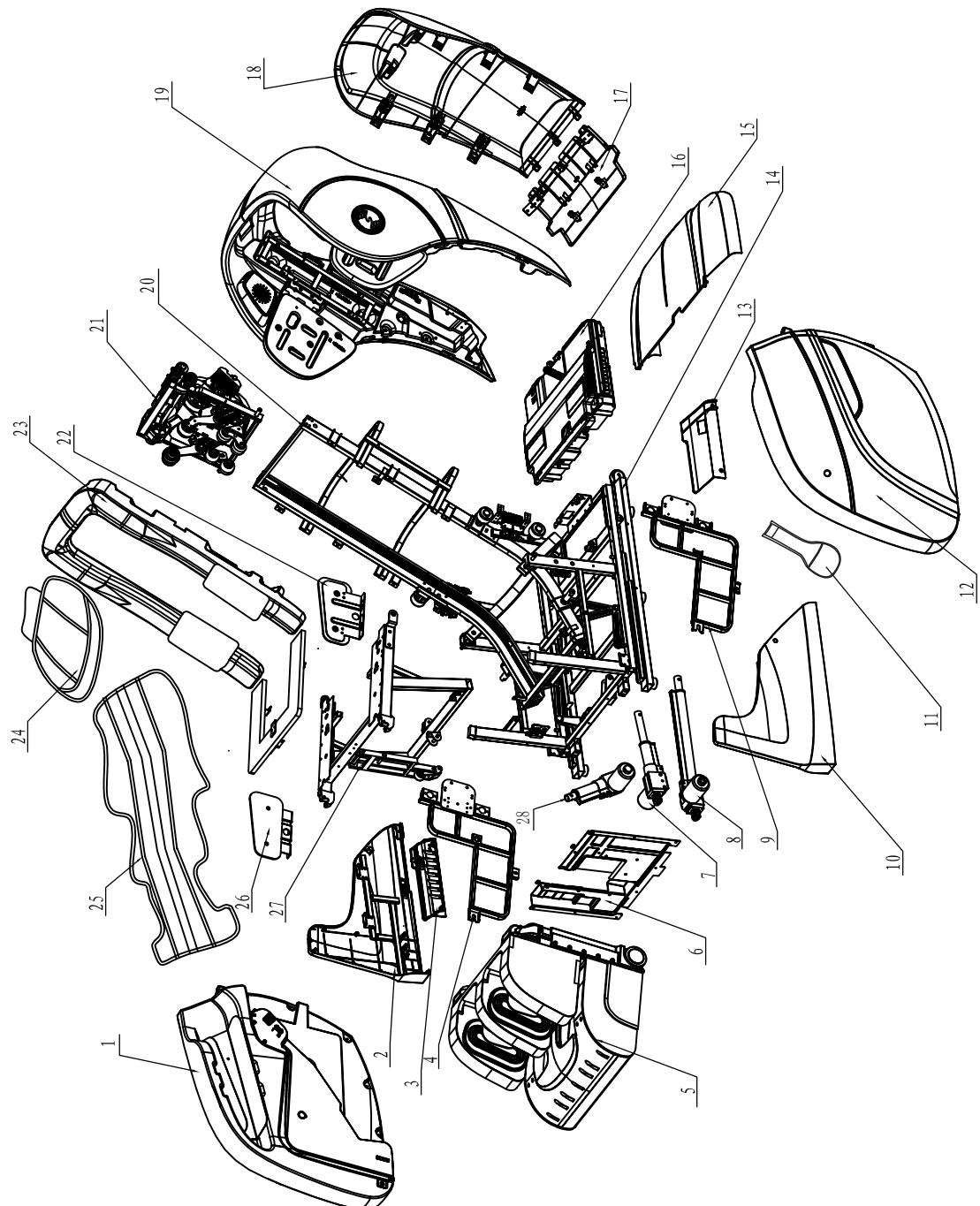

Mechanical components list

NO.	Name	Q'ty	NO.	Name	Q'ty
1	Right armrest assembly	1	22	Left side of seat airbag assembly	1
2	Right side fender	1	23	Backrest assembly	1
3	Right side of internal fender assembly	1	24	Pillow	1
4	Right armrest convert frame assembly	1	25	Back and seat cushion	1
5	Legrest and footrest assembly	1	26	Right side of seat airbag assembly	1
6	Front fender assembly	1	27	Seat frame surface assembly	1
7	Seat electric actuator	1	28	Footrest electric actuator	1
8	Sliding forward electric actuator	1			
9	Left armrest convert frame assembly	1			
10	Left side fender	1			
11	Remote control pocket	1			
12	Left amrest assembly	1			
13	Left side of internal fender assembly	1			
14	Seat bottom frame assembly	1			
15	Power box cover	1			
16	Power case assembly	1			
17	Back cover assembly A	1			
18	Back cover assembly B	1			
19	Shoulder device assembly	1			
20	Back frame assembly	1			
21	RK08 massage machine	1			

Installation (1)

Step 1. How to move the chair

- ① Take the chair out of the carton box including the chair body, back seat pad, pillow, power cord, remote control and control pocket.
(As the photo shows)



- ② Take out the left ,right armrest from the armrest carton box,Catching the two side of the armrest and take them out slowly. (The armrest is placed upside down in the box):put the armrest on the table or ground in flat(pay attention to prevent scratches), Catch the left and right two sides from the front of armrest and take it slowly.




Installation (2)

Step 1. How to move the chair

- ③ Take out the legrest device from the legrest carton box . Catching two sides of the legrest and take it out . Put it on the flat ground . (Pay attention to prevent scratches)

**Step 2. How to assemble the legrest**

- Install the legrest to the main body, ①Unscrew the bolt on the hanger of the legrest.
②Connect the legrest's air hose with the corresponding air nozzles of the main body tightly.
③Connect the two power cables' plug to the socket tightly , pay attention to the color and direction, black to black one, white to white one. ④ Install the rotation shaft bushings, raise the rotation shaft bracket on hangers, raise the legrest with hands, put rotation shaft into hangers on main body, and put rotation shaft bracket down. ⑤ Screw it back tight.
⑥Connect the legrest leather piece with seat by velcro .

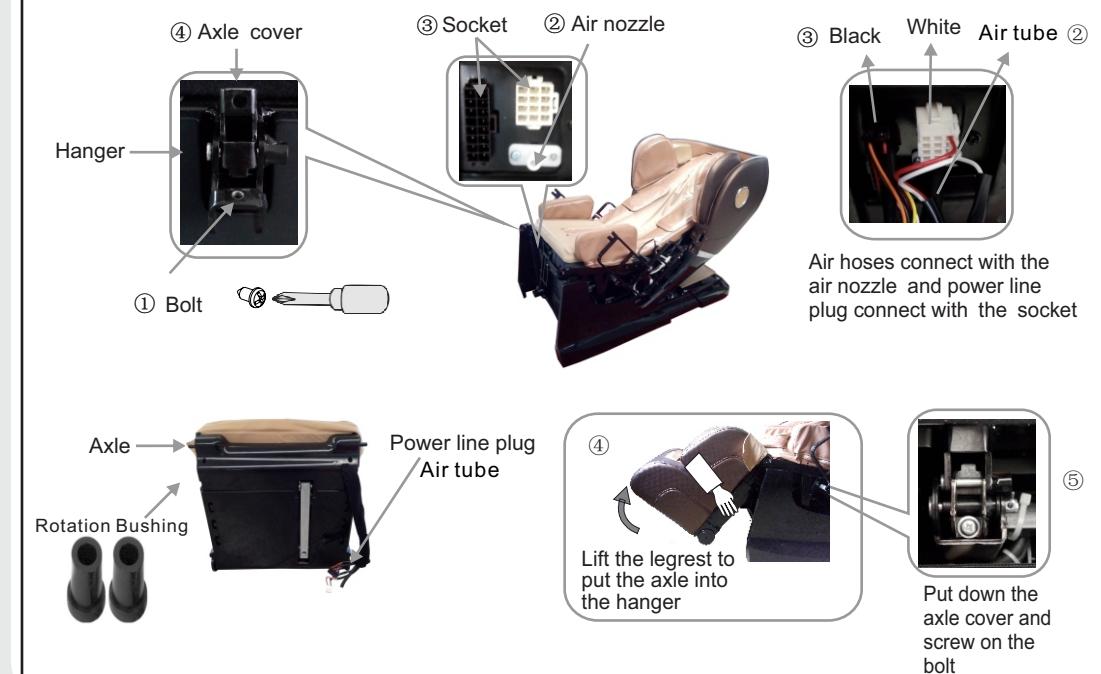

Installation (3)



Step 3. How to assemble the left and right armrest

Lift the front side of armrest and plug it into the rotating stand on the side of chair body. Connect the air hose of armrest with the air tubes on the side of chair body tightly. Black air tube should be connected with air hose marked with No. 1 and grey air tube with air hose marked with No. 2. Plug the power cord on the side of chair body into the socket of armrest in the right direction. ④Insert the rear side fixing plate into the fixing hole behind the chair body, and press the armrest downwards gently.

⑤Remove the rubber plug of armrest inner plate. Fix the screws inside the hole tightly, then put on the rubber plug tightly. ⑥Assemble the other armrest according to the above steps.

⚠ Don't squeeze air tube or power cord during the installation procedure.

◆ How to disassemble the armrest. 1.Adjust the back-rest to the lowest and restore the leg-rest before disassembly. 2.Remove the rubber plug inside the armrest and loosen the screw inside the hole. 3.Tilt the armrest a little upwards and disconnect the air tubes and power cord. 4.Continue to tilt the armrest upwards to dismantle it. 5.Insert the rubber plug into the round holes inside the armrest.

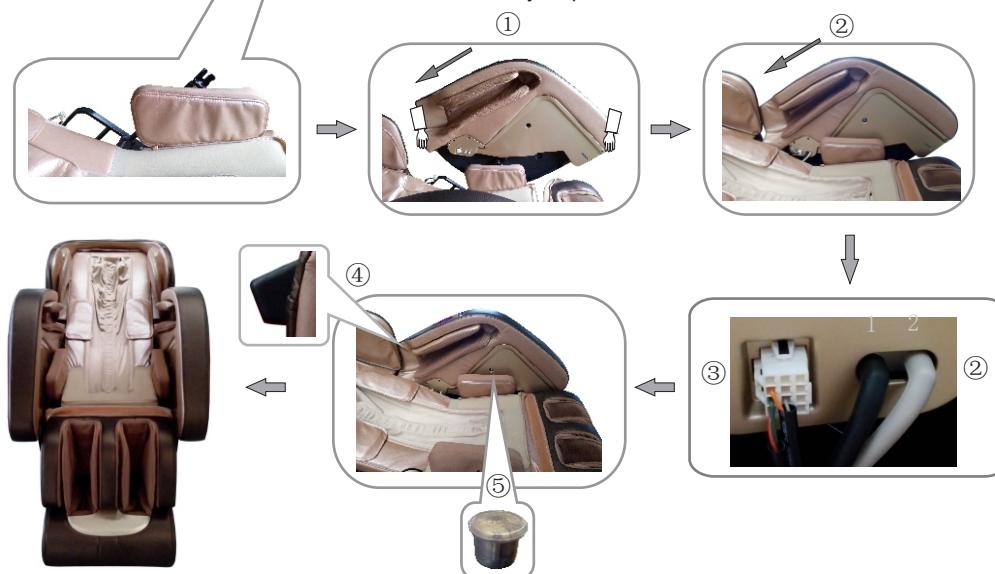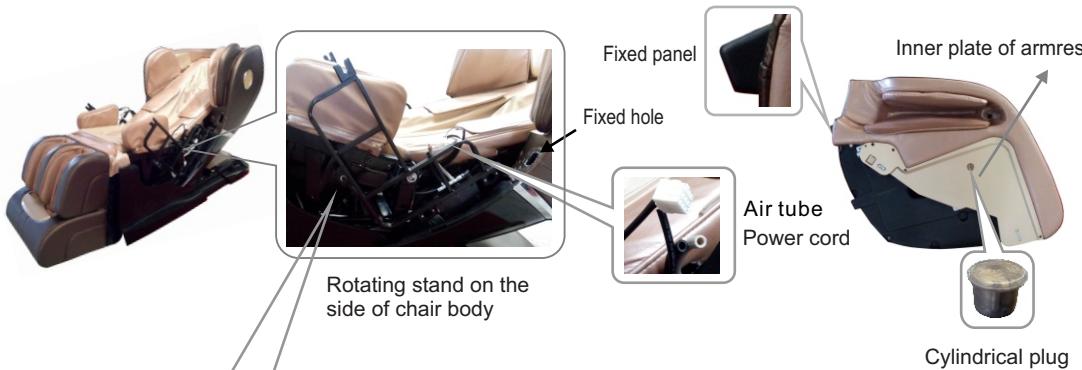


Installation (4)

Step 4. Fix the back seat pad, pillow and remote control pocket.

Put the back seat cushion, head pad onto the backrest and connect the two parts with zipper and velcro.


**★ Tips**

In order to get exact correct body test data for comfortable massage effect, when you sit on the chair, please make the bottom rim of the head pad flush with your ears before zippering the said two parts.

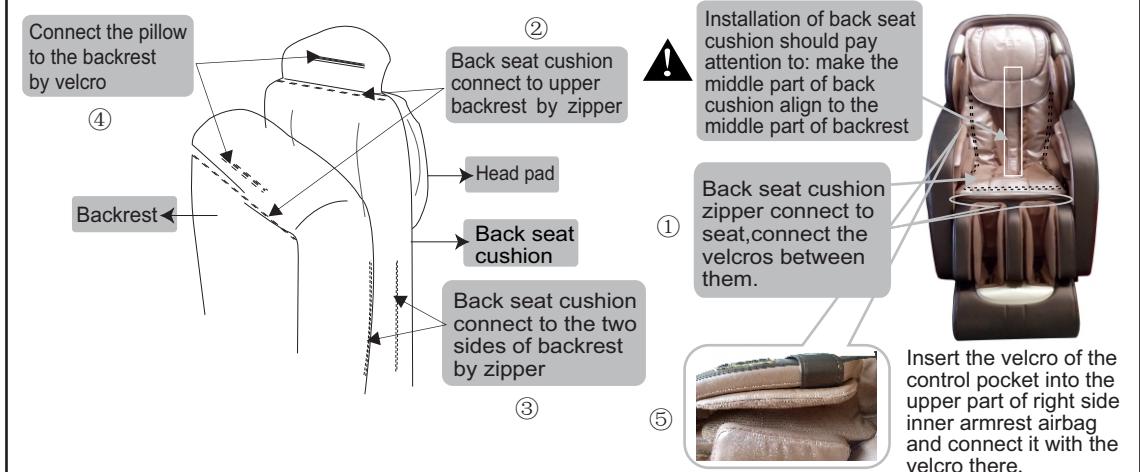**Step 5. Restore the chair.**

- ① Plug the cable line in the hole a on the power box as shown in the right picture.
- ② Insert the plug of cable line into the socket.
- ③ Turn on the switch behind the power box. (Switch on/off "O/I", press I down to turn it on.)
- ④ Press on/off button on the control twice to restore the chair. The chair will return to its initial state. Select the massage function by control or ipad. Enjoy the relaxing and pleasant massage.
- ⑤ In ON status, the USB port can be used for charging mobile.

⚠ Make sure the switch is off before you insert the plug of cable line into the socket.

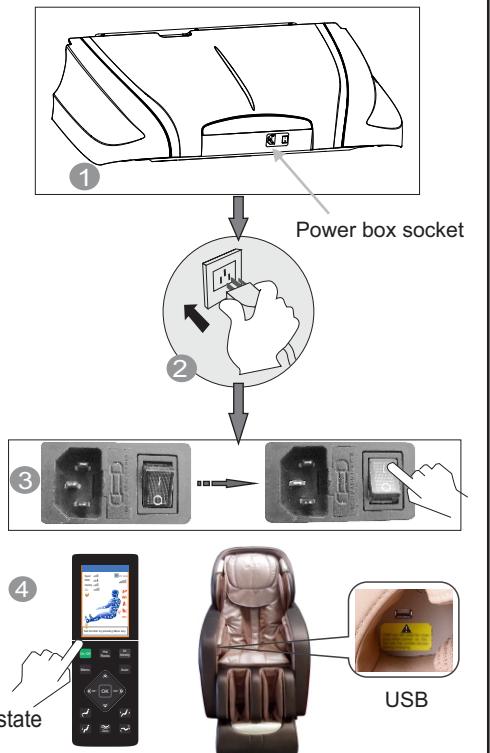

Environment**Caution**

Don't use the chair in wet places like near swimming pool, bath room to avoid leakage or electric shock.

Don't use under sunlight or near stove or other heating device to avoid upholstery deterioration

Please use chair in flat place to avoid chair falling or noise or other unexpected trouble.

Check places in front

Check the areas where the back and leg are extending. Make sure no people, pets or other things are around.

**Check power cable, plug**

- (1) Clean and dust from the plug. (using dry cloth).
- (2) Make sure there is no damage to the power cable.

**Check environment****Space size for using the chair**

- At least 4" from wall.



- At least 12" in front of the legrest.

Check surroundings

When the backrest or legrest reclinates or extends, please pay attention to:

- (1) Check any person or things in every direction.
- (2) Check the space size.

**How to move massage chair****Caution**


Don't move chair when it's working in case of accidents or damage

How to move

- ※ Bind the cables and put them in the seat to avoid damage when moving.
- ※ Do not drop the chair, as internal and external parts may get damaged.
- ※ Don't use the castor if the floor is wood or other easy damaged material.

**Normally**

- ① Lift up the legrest to make the gravity fall on the castor. Then push chair to destination and set it down.

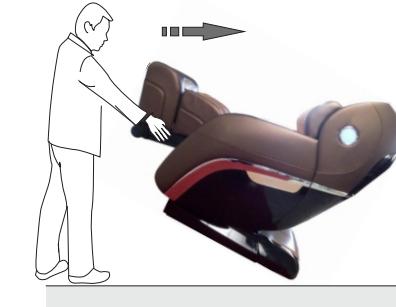

Please hold the footrest to move the chair as shown and be careful not to scratch leatherette.

**Adjustment before massage****Before adjustment**

- ★ Make sure the massage machine rollers stop in the upper part of backrest, then sit on the chair slowly to avoid wrong position of the massage head.

Body Scan

When doing body scan, if massage sensors can't detect the shoulder position, the chair default position is middle.



Right position: Your back clings to the chair back, your head clings to the pillow.

The remote control gives a "beep" sound after body scan is finished. When micro-adjustment is needed, just press up/down key on the remote control.

※ Make sure the body scan is finished correctly for the best massage.


Suggestion

- When you feel the intensity isn't strong enough, take off the pillow or cushion.

- Because the chair is heavy, it may mark floor. Please put a blanket or other protection down to insure the floor is not damaged.



Manual control

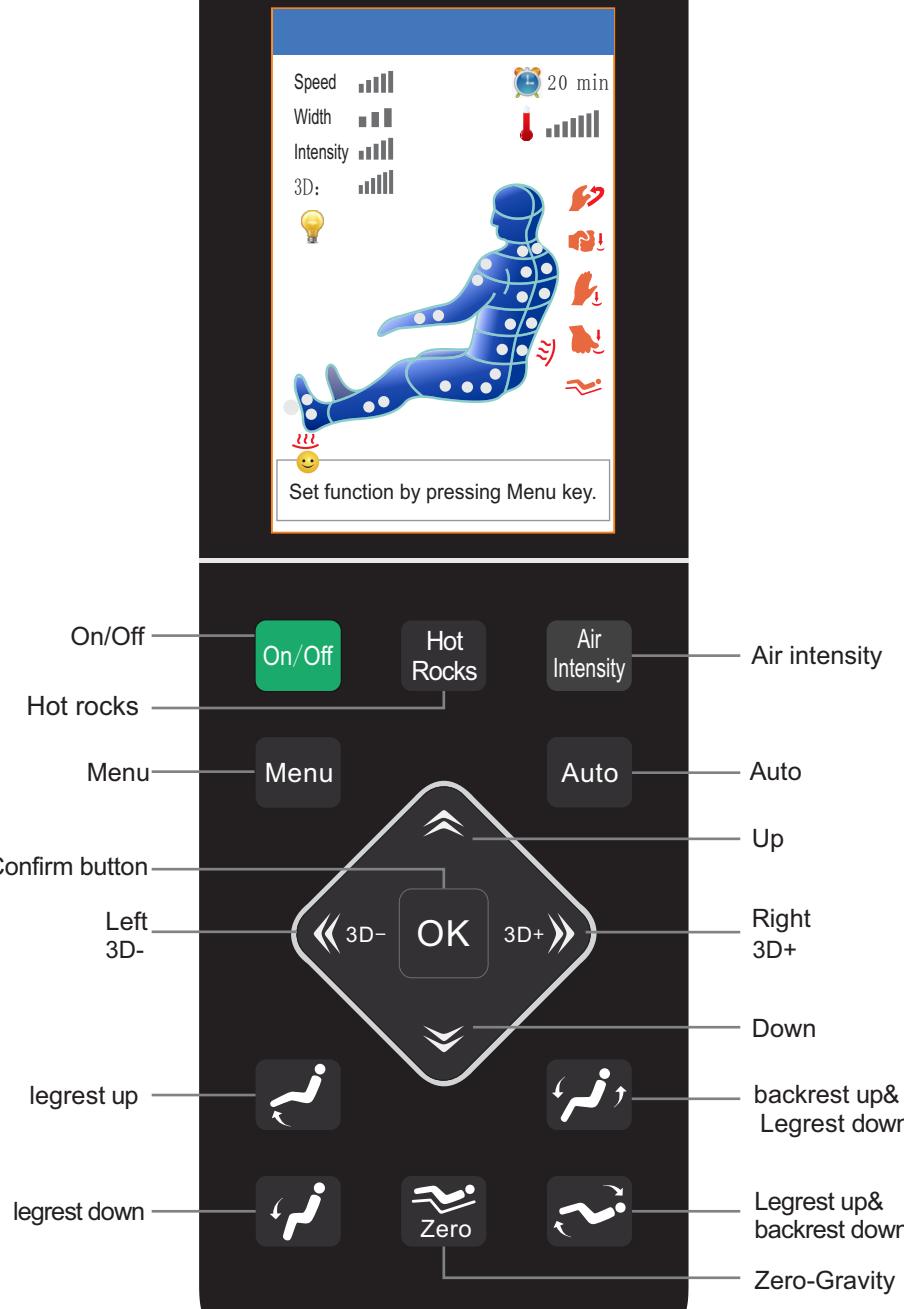

Massage operation

On/Off

Hot Rocks

Air Intensity

Menu

Auto




OK



1. Switch on, press on/off key to enter POWER ON status. Press Menu or other buttons to select any auto or manual program, massage machine will start the body scanning and then do the massage. In ON status, press on/off key to power off and get the chair back to upright position. Remote control will turn off in 3 to 5 seconds.
2. In ON status, press this key to enter HOT ROCKS program. Its default time is 5mins, click this button again, it will increase the time increment by 5mins, and the maximum will not exceed 30mins. In HOT ROCKS mode, massage machine heat is turned on and it will put on the max heat.
3. In ON status, press this key to adjust the air intensity either in auto mode or air massage program.
4. Press this button in Main interface to switch to menu interface, user can set different functions for the massage chair. Press it in Menu interface to return to main interface.
5. In ON status, press this key to activate different auto modes. Six common auto modes and six special auto modes for your choice, Recovery, Relief, Stretch, Music Sync, Renew, Calming, Refresh, Energy, Revive, Active, Relax and Healing. Relative screen display is shown on the controller.
6. There are four direction buttons, Up, Down, Left and Right. In Menu interface, it's a navigator button, user can press Up or Down to choose menu fast, and press Left and Right to enter previous or next menu. Besides, Up and Down button can be used to adjust the massage position in body detecting and point massage. In auto and manual mode, Left and Right button can be used to adjust 3D intensity.
7. In ON status and under Menu mode, press it to confirm menu choice and enter next menu.
8. Keep pressing this button, legrest goes up, release it, the legrest will stop moving. Controller screen will show position adjusting... when adjusting the position of the massage chair.
9. Keep pressing this button, backrest goes up, legrest goes down, release it, backrest and legrest stop moving. Controller screen will show position adjusting... when adjusting the position of the massage chair.
10. Keep pressing this button, legrest goes down, release it, the legrest will stop moving. Controller screen will show position adjusting... when adjusting the position of the massage chair.
11. Keep pressing this button, the backrest will go down, the legrest go up, release this button, the backrest and legrest will stop moving. When adjusting sitting position, it shows seat position adjusting... in controller screen.
12. In ON status, press this key, chair will go to zero-gravity position. Press it again to enter sleep mode. Press it once again to return to upright position. It will move like that in circles.

Massage operation

LCD display graphic description

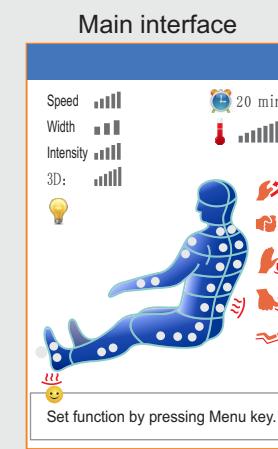


Press "on/off" button when power is on, a buzzer beeps once short and brittle. LCD screen displays boot screen, after 2-3 seconds it displays the main interface and enters standby mode. Main interface shows massage position, method, time, heating current massage status etc., all information will be updated timely when the chair is working.




- Massage position:**
- Kneading
 - Tapping
 - Knocking
 - Pressing
 - Zero gravity
 - Massage time left
- Massage function:**
- Massage speed (Level 1, 2, 3, 4, 5)
 - Massage width (narrow, mid, wide)
 - Air intensity (Level 1, 2, 3, 4, 5)
 - 3D (Level 1, 2, 3, 4, 5, 6)
 - Armrest LED light Status

Heat Status Display Temperature Level of Massage Machine Heat (Level 1, 2, 3, 4, 5, 6, 7)

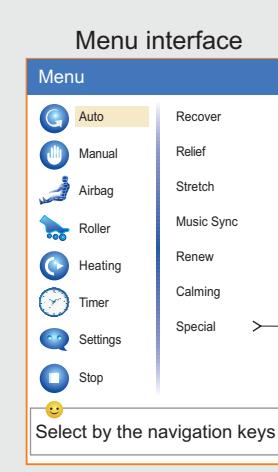

Press the "menu" button in the main interface, it will switch to the menu interface, user can set all massage function. Press the "menu" button in the menu interface will exit and return to the main menu interface.

Auto Mode: Press Menu to select auto mode by pressing OK. The special auto mode includes six special ones, Refresh, Energy, Revive, Active, Relax and Healing.

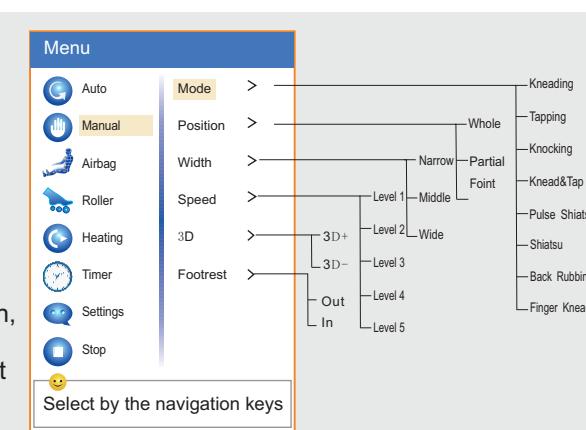


Manual Mode: Press Menu to select manual mode by pressing OK. It has six function settings including Mode, Position, Width, Speed, 3D and Footrest. Select one function, then press OK to enter next menu, select the desired massage mode and position, adjust the width and speed.

When massage chair in manual fixed point work modes, users can use up/down button to adjust massage position, information for "massage position adjusting..." will show at the bottom of the display screen.

Partial: Massage machine will center on the current position and move up and down in circles within certain limits.

- When select the massage mode, the default position is whole, just do only one time body scan.
- Partial and Point massage doesn't need body scan. In other words, if during the body scan or before the body scan completes, select Partial or Point will stop the body scan and directly enter the corresponding function, if you switch back to whole or select other mode at this time, body scan will be activated once again.


LCD display graphic description

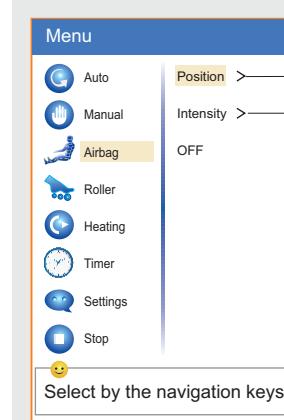

Airbag : After press MENU button, choose "air pressure" through navigation key. Press "OK" button , there are two kinds of air pressure related function: Position and Intensity ; Choose the corresponding function, then press "OK" to enter next menu and choose required function.


Roller: press "menu" button,select foot roller massage by the navigation key,press the "OK" button, Press "OK" button, choose rolling massage mode.

- During all the auto massage, the default setting for foot roller is on, you can press stop under Roller to turn it off.
- Under Airbag, when select Full Body or Legs& Feet, foot roller is on, press Stop under Roller to turn it off.
- When select Upper Body or Waist & Seat under Airbag, the foot roller is off, press Quick or Soft to turn it on.

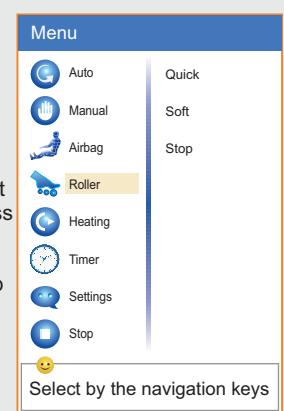

Heating: Press Menu and select Heat by navigation key. User can choose heat on or off. After heat is on, the heat symbol in the main interface will light up.

Timer: press "menu" key, through the navigation key to select menu item timing , can choose massage time: 30 minutes, 20 minutes, 10 minutes, the default time is 20 minutes, the chair will turn off automatically after time is up.

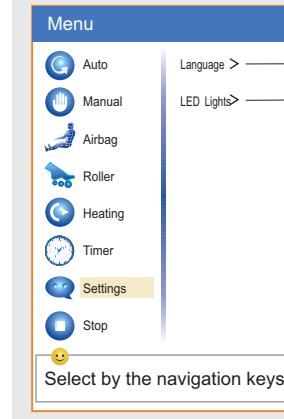

Settings: Press Menu, then set the language and LED light by the navigator buttons.


Stop: press "menu" key, select stop through the navigation key, it will stop massage, and go back to main interface.




After Power is on, the first time to use auto program,the massage chair will start body detecting, and "body scanning..." will be shown on screen display. After body scanning finishes, a buzzer will tick 5 times (once every second), meanwhile, the massage machine is located in shoulders position, if the user thinks the position is not proper, press Up and Down to adjust, when adjusting, "shoulder height adjusting..." is showing on screen display. Default working time is 20 minutes. (Massage chair will turn off and return after the massage is over)



ON/OFF

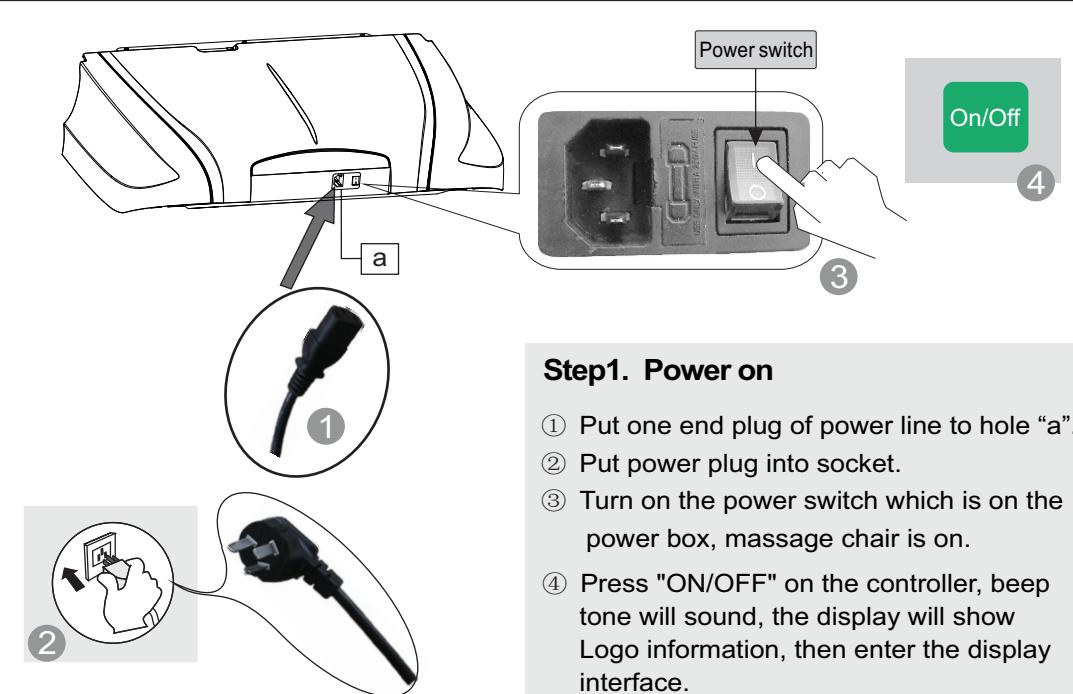**Step2. Massage**

- ① Operate according to controller operation instruction.
- ② When the massage stops, massage machine returns to the default position, meanwhile "Shutting down" is displayed on the controller, Controller turns off after 3-5 seconds.

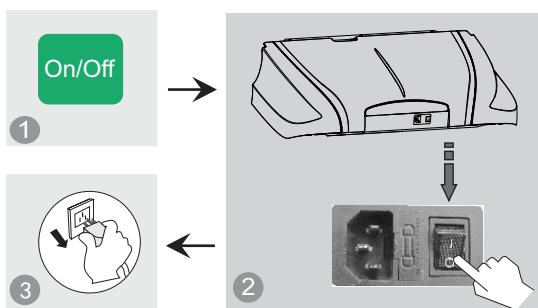

Warning	1. Before using, pull up the pillow, back cushion etc. to check if the PU material is broken, if there is damage, please stop using, unplug the chair, and call for repair. 2. Only plug in to proper socket.
----------------	--




1. Before sitting down, make sure there is no things pinched between armrest, legrest, backrest and seat of the chair.
2. Make sure the chair is in returned position before using.
3. Do not stand on it when using the massage chair.


Function Illustration

Legrest Up


Press this button on controller, legrest goes up. Release the button, legrest stops moving, the angle is locked in this position.

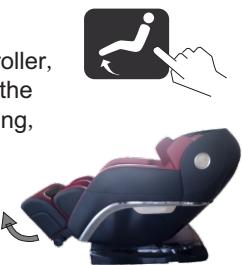**Legrest Down**

Press this button on controller, legrest goes down. Release the button, legrest stops moving, the angle is locked in this position.

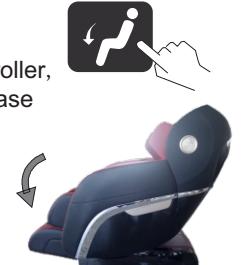**Legrest up and backrest down**

Press this button on controller, backrest go down and legrest go up. Release this button, backrest and legrest stop moving, the angle is locked in this position.

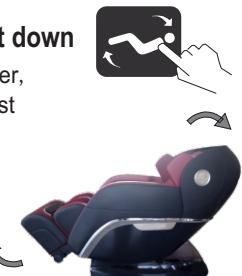**Backrest up and Legrest down**

Press this button on controller, backrest go up and legrest go down. Release this button, backrest and legrest stop moving, the angle is locked in this position.

**Zero gravity**

Press the remote control button ,the back goes down, the legrest goes up. The chair then goes into Zero gravity .Press the button again to reset the chair back from Zero gravity.

**Tips:** Zero-gravity

is based on the NASA Space capsule relaxation theory. Using the back ,seat, and legrest the chair can make the angle 128 degrees between Body and legrest. This puts the heart and knee at the same level, reducing the blood circulation pressure.,Not only that, but all the gravitational forces can be dispersed to the massage chair, reducing heart and spine pressure. This makes the Zero gravity position the most relaxing position possible.



Function illustration

Legrest extend and retract

Legrest can extend and retract within 8" automatically.

**Arm massage**

This Massage chair uses capsule type armrests, with built-in pulse air pressure massage functions. Put the arms between armrest's air bags. After filling with air, kneading massage begins on the arms and hands.

**Shoulder massage**

Multi-layer extrusion type air bags on both sides of the shoulder, realize stretch function. The 3D massage machine pushes your body while air pressure holds on to your shoulders, allowing for body stretch function.

**Warm foot therapy**

Unique to this chair, two-way rotary foot kneading technique, realizing the sole function of kneading foot, combining with foot heating. It can perform the best foot massage experience.



- 1. Please make sure there is no people, no pets or any things within range for backrest to avoid unexpected damage or property loss.
- 2. When adjusting the backrest, be careful of the gap between the backrest and armrest.
- 3. When the legrest is moving, please do not stand on or leave from the seat, do not stop the legrest from moving in case of causing damage to the chair.
- 4. Please refer to control instruction for more details.

**New Curved rail. Unique back and buttocks massage.**

Massage chair adopts 3D curved rail technology, combining with 3D massage mechanism patent, realizing multi-angle integrated massage for back and buttocks. It has all the benefits of traditional L shape rail function, as well as adding in the ability to lay down flat and fully stretch the body.

**Music therapy**

There are multiple speakers around your head. The blue-tooth wireless music player function integrates music into the massage techniques to realize a music sync relaxing therapy and effectively adjust the mood and ease fatigue.



After massage

Restoration of massage chair

- ① Make sure there is no people or pets or anything within the range of backrest and legrest.



- ② Make backrest and legrest restore to the normal position by controller.(or press ON/OFF key to make backrest and legrest restore to the normal position)




- ③ Adjust the back cushion and pillow to the original position.




- ④ After chair returns to the upright position, put the remote control back in the pocket.

- ⑤ Switch off the power.

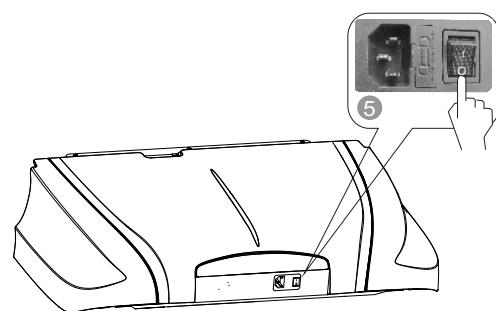

- ⑥ Unplug the electricity from the socket.

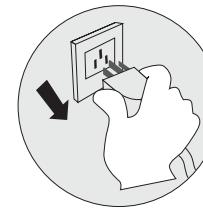

- ⑦ Clean the chair if necessary.



Attention!

- After using, press ON/OFF to end all the massage method and wait for restoration. Make sure to cut the electricity to avoid any unexpected hurt to children or pets. After using, turn off the power, unplug the electricity.

Synthetic leather**General cleaning**


Clean with tidy, soft and dry cloth.



(Chemical or medical cleaning products are forbidden)

**If Synthetic leather is dirty**

- ① Dip soft cloth into 3-5% neutral detergent solution, squeeze it then clean the stain on the surface gently.
- ② Wash cloth, squeeze it and wipe with detergent.

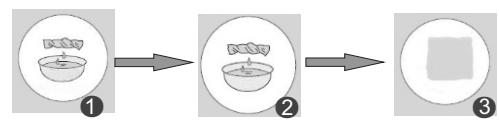

- ③ Clean it softly with dry cloth.
- ④ Allow to dry (do not blow dry)

Plastic parts

- ① Dip cloth into neutral detergent, squeeze it and clean.
- ② Wipe with water-dipped but squeezed cloth.
- ③ Allow to dry.

**Clean the controller**

- ① Wipe with soft cloth after squeezing.



- ② Allow to dry.

Propellant, diluent or alcohol are prohibited.

**Fabric material**

- ① Dip cloth into weak neutral detergent, squeeze it and clean the material surface.

- (Diluent, gas, alcohol are prohibited)

Brush the dirty parts with weak neutral detergent.
(avoid excessive brushing)

- ② Excessive cleaning may cause damage to cloth material.

- ③ Wipe with water-dipped but squeezed cloth.

- ④ Allow to dry.

**About cleaning**

- Before cleaning, unplug the electricity, don't touch plug with wet hands.
- Make sure hands are dry when plugging to the socket.
- Failure to follow these instructions may cause user to experience electric shock.

FAQ and Troubleshooting

Trouble	Possible reason	Troubleshooting
The massage chair won't work after power is on.	Check if the control is on or not.	Open control switch
	Not choose massage function	Choose massage function
	Check whether the plug and socket are in poor contact.	Make sure plug and socket are connected.
	Wire or plug is damaged.	Please ask manufacturer or appointed or authorized maintenance man for checking and replacement.
The fuse burned out.	Replace fuse of same specification.	
Something wrong with inner circuit	Please ask appointed or authorized maintenance man for checking and repairing.	
There is beep sound inside the massager while working.	It is sound of air pump, motor and mechanical parts.	Normal sound, no need to pay attention to it.
Different height of the left roller and right roller.	Massage rollers work on shifts.	Normal, no need to pay attention to it.
Sound suddenly becomes much louder while working.	Check if it is due to long time working.	Turn off the power to have the massager rest for over half an hour before next turn working.
	Inner parts worn out from long use.	Please ask appointed or authorized maintenance man for checking and repairing.
Massage function suddenly stops while working.	Sudden power-loss	Connect the electricity power again
	Setting time is up.	Turn off the power to have the massager rest for over half an hour before next turn working.
Neither the backrest nor the legrest can move.	Check if there are any obstacles; the chair has safety functions that will shut it down if there is an object blocking normal function	Remove the obstacles; turn off the power to have the massager rest for over half an hour before next turn working.
The chair can't return to its original position.	Excessive operation; setting time is up.	Please ask manufacturer or appointed or authorized maintenance man for checking and replacement.
Electric wire or plug hot.	It may have overload	Turn off the power to have the massager rest for over half an hour before next round working.
	Other abnormal case	Ask or appointed or authorized maintenance man for help.



Note: If the troubles are not caused by above reasons, please turn off the power switch and cut off the power, then contact the manufacturer or service agent to repair your chair.

Technical data

Name	INFINITY	Power	280W
Model	Presidential	NW/GW	Main body :253/302 lbs Armrest :32/43 lbs Legrest &footrest : 66/76 lbs
Input voltage	<input type="checkbox"/> AC120V, 60Hz <input type="checkbox"/> AC220V, 50/60Hz <input type="checkbox"/> AC240V, 50Hz	Packing size	Main body: 59.8×31×38 " Armrest : 38.5×15.5×28 " Legrest &footrest :23 × 20 × 24 " (L×W×H)
Safety design	Class I		

Others